

REVOLVER[®]

UNIVERSAL DRILL BIT GUIDE

The Revolver: Universal Drill Bit Guide[®] System revolutionizes conventional drill bit guidance by replacing multiple standard guides with one, 10-piece kit.

COMPONENTS:

- The Revolver: Universal Drill Bit Guide[®]
- Sterilization Tray
- 8 Tip attachments

For Small Diameter Bits

(1.0-2.5 mm) in 130mm, 160mm, 190mm, and 220mm lengths

For Large Diameter Bits

(2.6-5.0 mm) in 130mm, 160mm, 190mm, and 220mm lengths

FEATURES:

- Rotating chamber adjusts to multiple drill bit diameters from 1.0mm to 5.0mm
- Range of attachments extend the length of the guide up to 120mm for percutaneous applications
- Non-slip tip provides stability on the bone
- Comfortable, ergonomic handle
- Hand placement is safely away from the drilling point

BENEFITS:

- Cost savings by eliminating numerous conventional guides
- Ease of inventory management and sterilization
- Reduces surgical time in switching and locating various guides
- Little to no training required
- Universal use with most standard drill bits



60mm x 2.5mm



90mm x 2.5mm



120mm x 2.5mm



30mm x 2.5mm



30mm x 5mm



120mm x 5mm



90mm x 5mm



60mm x 5mm


REVOLVER[®]

UNIVERSAL DRILL BIT GUIDE

The Revolver Universal Drill Bit Guide[®] can be used with all orthopedic surgical drills. It protects soft tissue and grips to the bone to promote stability during drilling.

"With a turn of the chamber, I can switch drill bit guide size with ease. It saves me time and saves the hospital money."

- Orthopedic Surgeon, Philadelphia, PA

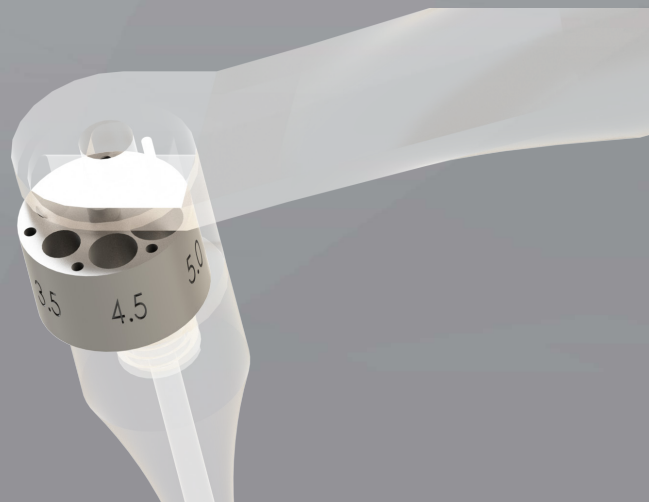

For more information, or to order your Revolver Universal Drill Bit Guide[®] directly, call 307.315.6403 or visit www.mcginleyorthopedics.com.

ALSO FROM MCGINLEY ORTHOPEDICS

IntelliSense Drill Technology[®]

- **Real Time Depth Measurement** - Eliminates the need for a depth gauge, displays the depth of holes instantaneously
- **Auto-Stop Technology** - Drill bit rotation stops when it passes through the far bony cortex, auto-stops at a designated location
- **Edge Detection** - Alerts surgeon when passing through any bony cortex, signals dense bone/uncortical breach
- **In-Sight Surgical Lighting™** - Directly illuminates the surgical field using integrated LED lights



Learn more at www.intellisensedrill.com/drill. Call to order today! 307.315.6403

www.mcginleyorthopedics.com